

ElvaX S Lab


ElvaX S Lab is designed especially for analysis of petrochemicals for sulfur content in wide concentration range.



Due to low detection limits (down to ppm level) and minimal required sample preparation ElvaX S Lab is an ideal tool for analysis of sulfur in petrochemicals.

Many international standards describe analysis method for sulfur in fuel and oil. Sulfur forms damaging sulfates in vehicle exhaust and pollute the atmosphere by sulfur oxides.

Furthermore sulfur reduces fuel stability and effects on engine performance. Maximum permissible level of sulfur in fuel decreases from year to year.

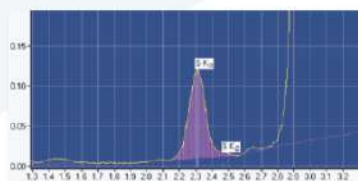

Lowest detection limit

Unique configuration of x-ray tube and geometry designed for best sensitivity to sulfur.
Limit of detection is only 2 ppm!



Fast analysis

Sample preparation is not required for analysis of petrochemicals. Result will be obtained after one minute.



One instrument for any sample

Analysis of gasoline, oils, lubricating, diesel, crude oil and any other petrochemicals using one instrument. No additional calibration is needed.



Leakage protection

The movable analyzer stage protects expensive components from the leakage of the investigated petrochemicals.



ASTM D4294 and ISO 8754

2003/17/EC (or EURO VI) directive requires a maximum amount of sulfur in fuel of 10 ppm. ElvaX S fully meets the requirements of ASTM D4294, ISO 8754, ISO 20847, IP 336, ASTM D6445.

X-ray tube
Anode: Ag
Wattage: 4W
Voltage: 40kV

X-ray Detector
Type: Fast SDD
Area: 25 mm²
Energy resolution: 140 eV at Mn Ka, 90 eV at Al Ka
Count rate: 500 000 cps

Electronics
Signal processor: Digital pulse processor based on 80 MHz DSP
ADC: 4096 channels

Common
Analytical Chamber Dimensions: 185 mm × 212 mm × 90 mm
Weight: 7 kg
Power supply: 90 – 240 V, 50/60 Hz; optional battery
Power consumption: 20 W

Embedded PC
Operating system: Windows © Compact
Analysis algorithms: Regression algorithm, Manual spectra comparison
Data transfer: 2×USB, MicroSD, miniUSB
Screen: 5", 800×480

Taqnia General Trading LLC

P.O. Box: 301, Office No. 3,
Al Rashida Building, Al Muwaihat,
Ajman-UAE.

Phone: +971 6 746 3141
+971 56 995 8998

Email: info@taqniatrading.com
sridhar@taqniatrading.com

www.taqniatrading.com